

**L&L SKIN**

**KIRA 24K GOLD RF FACIAL DEVICE**

**USER MANUAL & WARRANTY INFORMATION**



### **Introduction**

Thank you for purchasing your new RF Beauty instrument. The new radio frequency beauty instrument uses advanced radio frequency technology to target the epidermis, dermis, and muscle layers for multi-dimensional anti-aging, faster heating, larger effect area, and richer effect levels! Upgrade the 3-band frequency conversion radio frequency, activate collagen at multiple levels, make the skin plump and elastic, lift the wrinkles, and improve the four major aging problems of loose, sagging, wrinkled and yellow skin.

At the same time, it adopts dual ion penetration and purification technology, combines high-end ion technology and EMS technology, and acts simultaneously with micro-vibration massage to deeply purify the skin, help the skin to better absorb the nutrients of skin care products, and lift and firm the skin. Get better skin care results from the inside out and keep your skin looking healthy, supple, and radiant at all times.

### **Important**

Read this user manual carefully before you use the instrument and keep it for future reference.

### **Danger**

- Keep the charger and USB cable dry. Never immerse the into water.
- Do not charge this instrument in wet surroundings (e.g. near a filled bath, a running shower or near a filled swimming pool).
- Never throw the instrument into fire.

### **Warning**

- Only use the instrument in combination with the USB cable supplied when charging.
- Do not use the instrument, the attachment or the USB cable if it is damaged or broken to avoid injuries.

- If the USB cable is damaged, always have it replaced with one of the original type in order to avoid a hazard.
- Fully charge the instrument at least every 2 months to sustain battery lifetime.
- This instrument is not intended for use by persons younger than 12 years and persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the instrument by a person responsible for their safety.
- Women who are pregnant or lactating are not recommended to use this instrument.

### **Caution**

- Do not use the instrument for any other purpose than described in this manual.
- Charge, use and store the instrument at a temperature between 5°C/40°F and 35°C/95°F.
- Do not use the instrument if you are taking steroid based medication.
- Do not use this instrument if your skin is still recovering from a cosmetic procedure, such as a chemical or mechanical peeling, laser resurfacing, fractional laser treatment, or mesotherapy treatment or any other treatment with dermal fillers in the past two weeks..


### **General**

- Regular cleaning of the instrument ensures optimal results and a longer life of the instrument.
- Do not use when sleeping or drowsy.

### **Contraindications and other warnings**

- Do not apply this instrument directly to your ear as this may cause temporary hearing loss.
- Do not use this instrument if you have any ear problem (Swimmer's ear, ear infection, any anomaly).
- Do not use the instrument on chapped skin, open wounds, or if you suffer from skin diseases or skin irritations, such as severe acne, sunburn, skin infection, etc.
- Do not use this instrument if you have a pacemaker, implanted instrument, or have a metal prosthesis.
- Do not use the instrument on the sensitive skin area around the eyes, as it is thin and delicate.
- Do not use the instrument on any mucosal area (e.g. inside nose, mouth, genital area) as it may cause local irritation.


## Product Mix

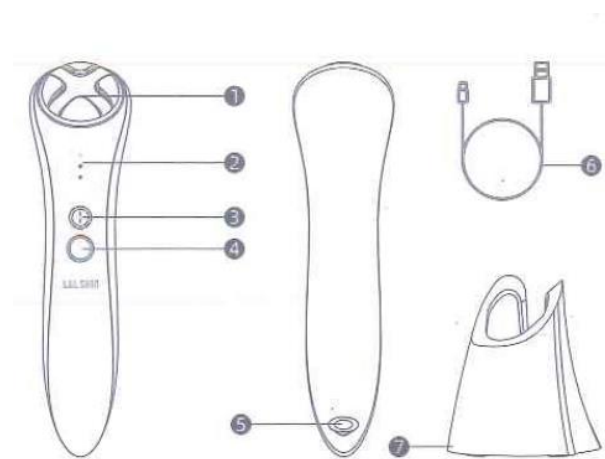

1. Massage Probe
2. Gear Indicator/Charging indicator
3. Gear Switch Button
4. Power Button, Mode Switch, Mode Indicator
5. Charging Socket
6. USB Cable
7. Base

## Function Guide

1. The instrument consists of four modes: skin cleansing mode, nutrition mode, lighten wrinkles mode, and firm skin mode.
2. After long pressing the power button to turn on the instrument, the four modes can be switched by the mode switch button.
3. After the instrument is turned on, it can be switched by the gear switch button.
4. Press and hold the power button to power off.

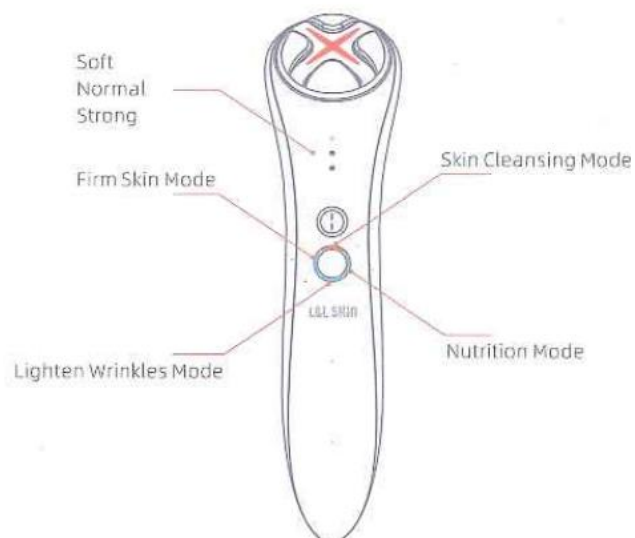

\*Auto shut off function: the instrument will shut off automatically after 4 minutes of continuous use.

## Function Guide

### Skin Cleansing Mode

1. In this mode, the continuous frequency conversion vibration and positive ions generated by the instrument can purify the deep pores, shrink pores, lighten spots and acne, solve dull skin, and make skin more crystal clear. The specific wavelength of blue light can fight acne and purify the skin, shrink pores, repair sensitive skin, balance oil.

2. How to use

**Step 1:** Press and hold the power button to turn it on, and then select option to activate the skin cleaning mode.

**Step 2:** After applying the nourishing cream, touch the facial skin with the massage probe as shown in the picture and massage in circular motions.



### Nutrition Mode

1. In this mode, the continuous frequency conversion vibration and negative ions generated by the instrument can promote the skin to absorb nutrients, such as vitamin C and minerals. Therefore, when using this mode, the effect of using the instrument with skin care products is more obvious, such as whitening lotion, containing vitamin C, vitamin A or vitamin E skin care products.

2. How to use

**Step 1:** Press and hold the power button to turn it on, and then select option to start the nutrition mode.

**Step 2:** After applying the nutritional cream or mask, touch the facial skin with the massage probe as shown in the picture and massage in circular motions.



### Lighten Wrinkles Mode

1. In this mode, the radio frequency beauty instrument emits radio waves into the subcutaneous tissue, so that the natural resistance movement of the subcutaneous generates heat energy. It uses the principle that the collagen fibers in the dermis will shrink immediately at 55-70°C, which has the effect of tightening the skin and prompting collagen regeneration. Two effects, while the collagen produces immediate contraction, it will also stimulate the dermis to secrete more new collagen fibers to regenerate, thereby supporting the skin again, increasing the thickness and density of the dermis, and filling in wrinkles, improve relaxation.



The LED red light power adopts a specific wavelength, which can improve the microcirculation of the blood system and lymphatic system of the skin, stimulate the activity of mitochondria in the cell, and the specific wavelength of the LED red light can convert light energy into cell energy and neutralize the free genes generated in the cells, stimulate fibroblasts to produce collagen, activate the skin to speed up blood circulation, and can improve wrinkles and emollients from the inside of the skin.

The continuous frequency conversion vibration generated by the instrument can massage the skin, dilute fine lines, firm the skin, and effectively prevent wrinkles and sagging skin.

**Warning:** The radio waves generated by the radio frequency will continue to heat the dermis layer of the skin. Please continue to move the massage while the massage probe is attached to the skin. DO not stay in the same area of the skin to avoid burns.

## 2. How to use

**Step 1:** Press and hold the power button to turn it on, and then select option to activate the lighten wrinkles mode.

**Step 2:** After applying the nutritional cream, touch the skin with the massage probe as shown in the picture and move it upwards to massage.



## Firm Skin Mode

1. In this mode, the instrument sends out EMS micro-current to rejuvenate facial muscles and motor nerves, and promotes more effective contraction of facial skin, thereby promoting local blood circulation, activating cell, reducing skin puffiness, and lifting and firming the skin.

The LED red light power adopts a specific wavelength, which can improve the microcirculation of the blood system and lymphatic system of the skin, stimulate the activity of mitochondria in the cell, and the specific wavelength of the LED red light can convert light energy into cell energy and neutralize the free genes generated in the cell, stimulate fibroblast to produce collagen, activate the skin to speed up blood circulation, and can improve wrinkles and emollients from the inside of the skin.

The continuous frequency conversion vibration generated by the instrument can massage the skin, dilute fine lines, firm the skin, and effectively prevent wrinkles and sagging skin.

## 2. How to use

**Step 1:** Press and hold the button to turn it on, and then select option to activate the firm skin mode.

**Step 2:** After applying the nutritional cream, touch the skin with the massage probe as shown in the picture and move it upwards to massage.



## How to charging

1. Connect the computer to the charge with a USB charging cable. Turn on the computer, connect the USB charging cable to the computer charging socket, and connect the Type-C end of the charging cable to the charging socket of the instrument to start charging.
2. Charge with adapter: Connect the USB charging cable to the adapter socket, and connect the Type-C end of the charging cable to the charging socket of the instrument to start charging. Adapter specifications: 5V-1A
3. The indicator flashes when charging, and the indicator stays on when fully charged.




**CAUTION:** Before charging, make sure the charging socket is dry and free of water and skin care products. DO not charge the instrument where it may be submerged in water. Please do not use the instrument while charging, and the charging time should not exceed 24 hours. Fully charged the instrument at least every 2 months to preserve battery life. Please only use the charging accessories provided with the instrument.

## Cleaning the instrument

The following maintenance instructions are important to ensure the proper functioning of your instrument. If these instructions are not followed, your device may provide less than ideal treatment or the device may stop working altogether.

Before start cleaning, make sure your instrument is powered off.

- When cleaning the instrument, only use a dry cloth.
- DO not use steel wool, abrasive cleansers, or corrosive liquids such as alcohol or acetone to clean the instrument.
- Thoroughly clean the instrument and probe after each use to ensure performance.
- Do not rinse your instrument under the faucet or immerse it in water.
- Make sure all parts are dry before reconnecting to the charging cable.

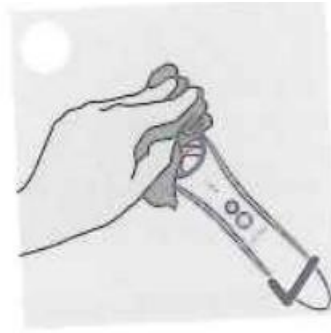

BATTERY: Li-ion 900 mAh 3.7 V

CHARGING: About 3 hours

Theoretical usage time for a full charge:

1. Lighten mode-about 40 minutes (depending on the skin resistance value will vary)
2. Other modes – about 300 minutes

STANDBY: About 60 days