


**L&L Skin**

**MIGA Portable Oral Irrigator**

**User Guide**

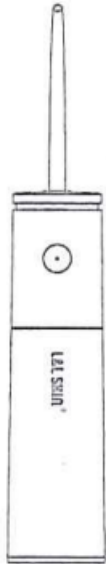

Thank you for choosing L&L Skin Portable Oral Irrigator.

Healthy teeth, happy life!

**Important safety measures**

When using this product, especially with children beside you, basic safety precautions must be taken.

**PLEASE READ ALL INSTRUCTIONS BEFORE USING**

**Dangerous:**

To reduce the risk of electric shock:

- Do not handle the charger with wet hands.
- Keep the product and charger away from water or other liquid.
- Do not reach for a product that has fallen into water. Unplug immediately
- Please check the product's intactness before use.

**Warning:**

- To reduce the risk of burn, electric shock, fire or other injuries.
- If the product has been opened or disassembled for any reason, the warranty will be invalid.
- Do not plug this device into an adapter that does not meet the specified voltage.
- Never operate this product under these circumstances:

It has a damaged cord or plug;  
It is not working properly;  
It has been damaged, or dropped into water.

Return the product to a service center for examination and repair.

- Do not direct the water flow under the tongue, into your ears, nose, or other delicate organs. The pressure produced by the product may do serious harm to these organs.


Please follow the operation steps when using this product.

- Magnetic fields: Pacemakers and similar implanted medical devices may be affected by magnetic fields.
- A slight bleeding is normal when using oral irrigator for first time. If there is heavy bleeding, stop using it and consult your dentist.
- Unplug this product before filling. Fill the reservoir with water only unless your dentist recommend other material. Never drop or insert any object into any opening or hose.
- Do not put the charger and wire on hot surface.
- Do not use this product around the place where aerosol (spray) products are being used or where oxygen is being administrated.
- Do not use iodine, slat water (salt solution or water-insoluble concentrated essential oil in this product. The usage of such substances will reduce product performance and shorten product life. If these items are used, it may invalidate the warranty.
- Do not use it while wearing oral accessories, and please remove them before use.
- Do not use it if you have open wounds on your tongue or mouth.
- If your doctor suggest that you receive antibiotics before surgery, you should consult your dentist before using this instrument or any other oral hygiene auxiliary chemical.
- Close supervision is necessary when this product is used by, on, or near children or invalids.
- Be sure to prevent children from playing with this electrical product.
- There are no user-maintainable parts in this appliance, and no maintenance is required.
- Only fill the water tank with warm water below 30°C or other solutions recommended by your dentist.
- Please use this product in accordance with the user guide or recommendations by your dentist.
- Never block air openings of the product or place it on a soft surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair, and etcetera.


Please keep the user guide.

If there are typographical errors or misunderstandings in the content of this manual , our company reserves the right of final explanation; in case of product upgrades or date updates, it will be incorporated into the new version of manual without further notice; if the appearance and colour of the product vary, the actual package product shall prevail.

## Product description

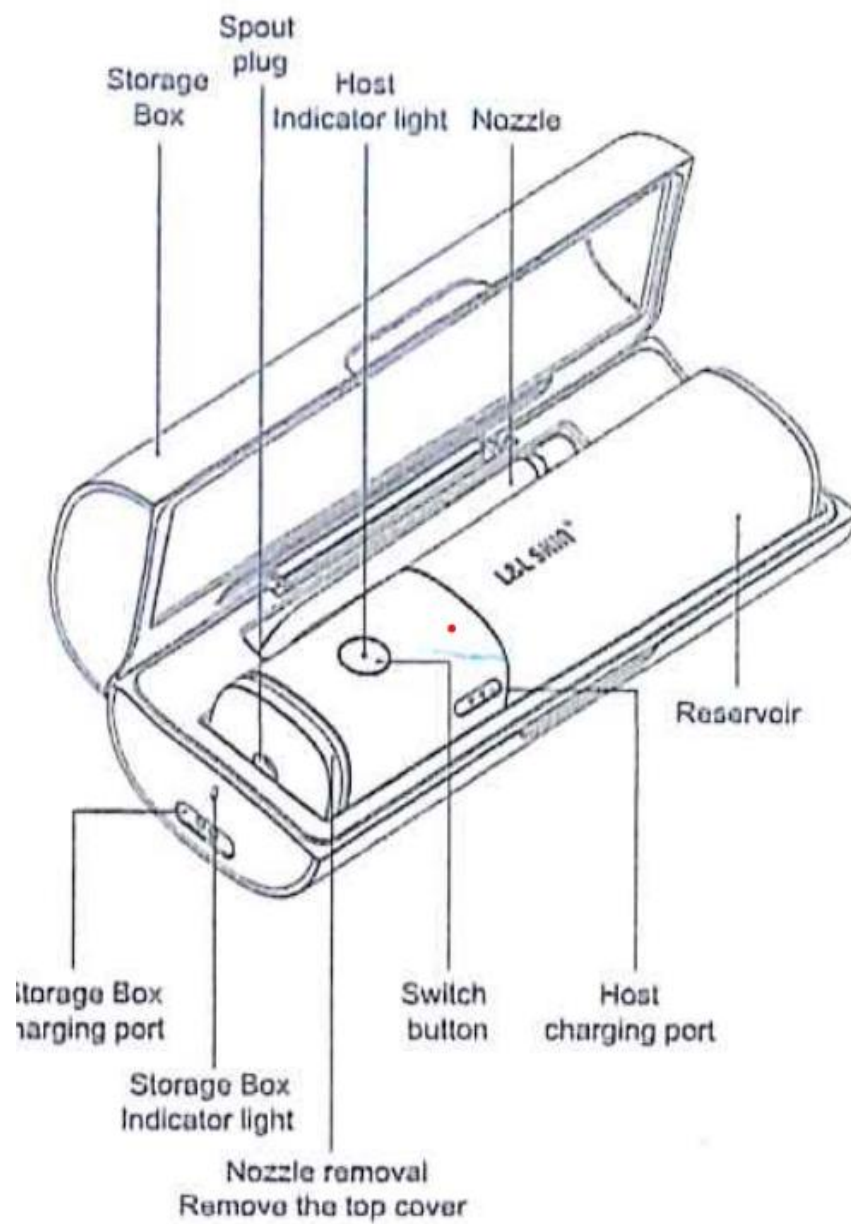

## Product accessories and specifications

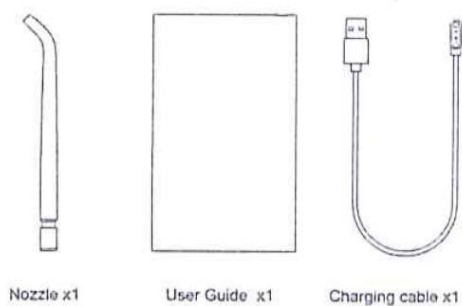

Replaceable nozzle/accessories can be purchased from the vendor.


### Product specification

Product name		L&L Skin Portable Oral Irrigator
Oral irrigator	Product number	MIGA
	Rated working voltage	3.7V
	Rated power	5W
	Waterproof level	IPX7
	Charging time	Requires more than 2 hours to fully charge
	Input parameters	5V 1A
Country of origin		P.R.C.

### Operation steps


#### Activation

**Very important!** Please fully charge the device before using it for the first time, connect the magnetic power cord to the charging port of the host and turn on the device. The product cannot be used during charging.

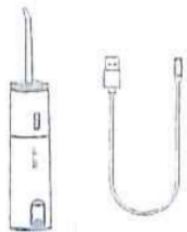


As shown in picture 1, open the reservoir cover and then pour water into it:

Note: Please make sure that the device is turned off when filling the reservoir.

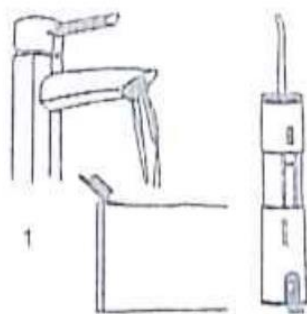

### Nozzle installation and removal

Insert the nozzle in the middle hole of the host until you hear the click sound, and the installation is successful. If you replace or remove the nozzle, please press down the nozzle to remove the top cover, and then take the nozzle out.

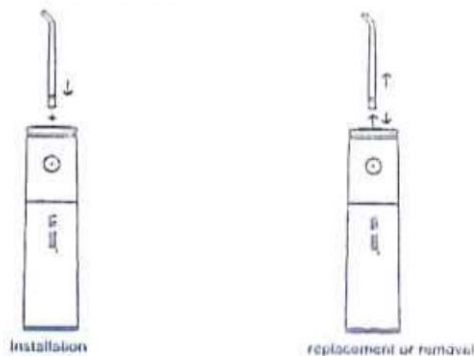

### Operation tips


The normal toothbrush can't clean the residue hidden deeply in your mouth. Traditional dental floss can only reach limited areas under the gums. But L&L Skin Portable Oral Irrigator can quickly and effectively clean food residue deep in the gums and between the teeth, keeping your mouth neat.

### Operation steps

- Press the switch button, the device starts to work, press it again, the device stops working. The device would power off automatically if it has not been operated in a minute.
- When charging, the magnetic power cord can be directly connected to the host for charging, or the host can be charged in a storage box.
- In the process of charging, the LD light flashes while charging and is on when fully charged.
- When using this product, the LED light keeps on with more than 20% power left and starts flashing or breathing below 20%.
- When the power is severely low, the LED will flash 3 times when you try to turn it on. It won't start.
- The charger is suitable for: AC 100-240VAC, 50/60 Hz, output 5V 1A.
- The device cannot be used while charging
- When in use, keep it upright to ensure the equipment is in the best working condition; by rotating the nozzle to adjust the water flow.
- Gum bleeding may occur the first time you use the oral irrigator, because everyone withstand different pressures in their gums. This is normal, so please feel safe to use it!
- When using the product for the first time, please start at the lowest gear and gradually increase the pressure as your oral cavity adapts.


If the water tank cover falls off:

1. Align the gasket of the water tank cover with the water tank mouth;
2. Press it in until you hear a "click" sound.

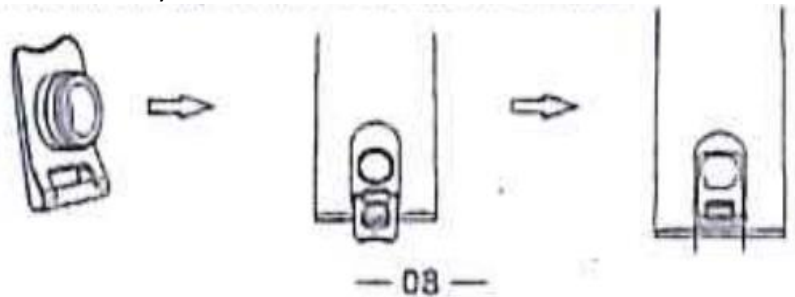


### Tips

Lean over the sink, put the nozzle in the mouth and aim at the gums and teeth, press the switch button to turn it on.

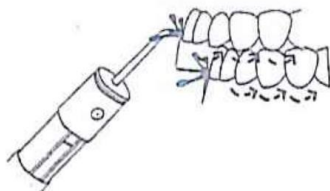


Spray the water at a 90-degree angle against the gums, and close your mouth slightly to prevent the water from shooting sideways. To archive the best results, please start cleaning from the molars (your back teeth) and moving toward the front teeth. Slide nozzle along the gum line and stay for a while between each tooth until the gums and the space between the teeth are cleaned.

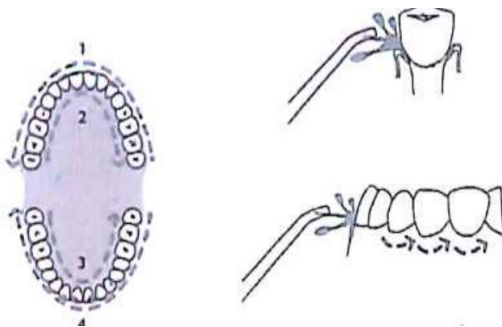


### Disassembly example

#### When finished using:

Turn off the device, open the back cover of the water tank with the injection hole facing down as shown in the picture, clear the left liquid in the water tank, turn on the device, and drain the left liquid form the device until no more water is discharged, then turn off the device.

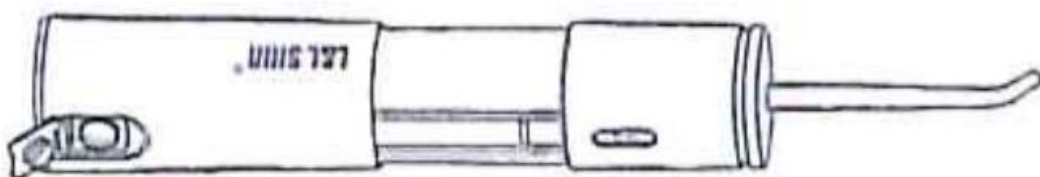

As shown in the picture, hold the product body with one hand and the water tank with the other hand while pulling the water tank down.

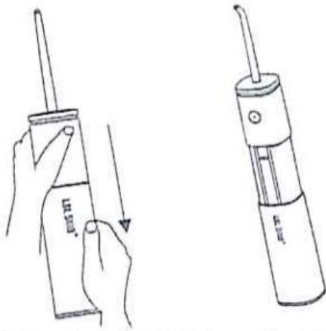

**Note:** There may be moisture inside the machine before its first use (in the reservoir and suction pipe), which is residue of the distilled water used for product tasting.

## **Battery handling and maintenance**

### **Battery maintenance**

- The device must be fully charged before first usage, and subsequent charges take about 2 hours;
- Plug the USB cable into the charging socket and connect it to a power source;
- If the device will be stored for a long time, please fully charge the device before storing;
- In order to maximize the service life of the battery, please charge it immediately when the power is low.

### **Battery handling**

This device contains a lithium battery, which should be recycled at the end of its life

Please recycle the battery at a local recycling center

The battery cannot be replaced.

### **Removal of the built-in battery**

- Use a screwdriver or other tools to pry and remove the nozzle to remove the top cover assembly, and uninstall the charging port.
- Remove the battery for recycling.
- Be careful not to short-circuit the positive and negative ends of the removed battery, and insulate the ends with tape.

## FAQ

Problems	Reasons	Solution
Leakage between the nozzle and the body of the device	The seal ring in the nozzle holder is damaged, the nozzle is not installed properly	Contact the vendor for repair, install the nozzle properly
Leakage in the water pump inlet	The seal ring is damaged, or there are impurities in the inlet side	Contact the vendor for repair
Leaking water tank	Damaged seal	Replace the seal ring, contact the vendor for repair
No enough pressure	Low voltage in use	Use after fully charged
Device cannot start	Low voltage in use	Use after fully charged or contact the vendor for repair
The product does not spray water	First use or not used for a long time	After adding a proper amount of water to the nozzle installation port, reinstall the nozzle, and turn on the power switch.

If you still can't solve the problem according to the above methods, please consult the vendor shop or go to a designated after-sales repair shop for repair.